


# Uretero-Renoscopes (URS)

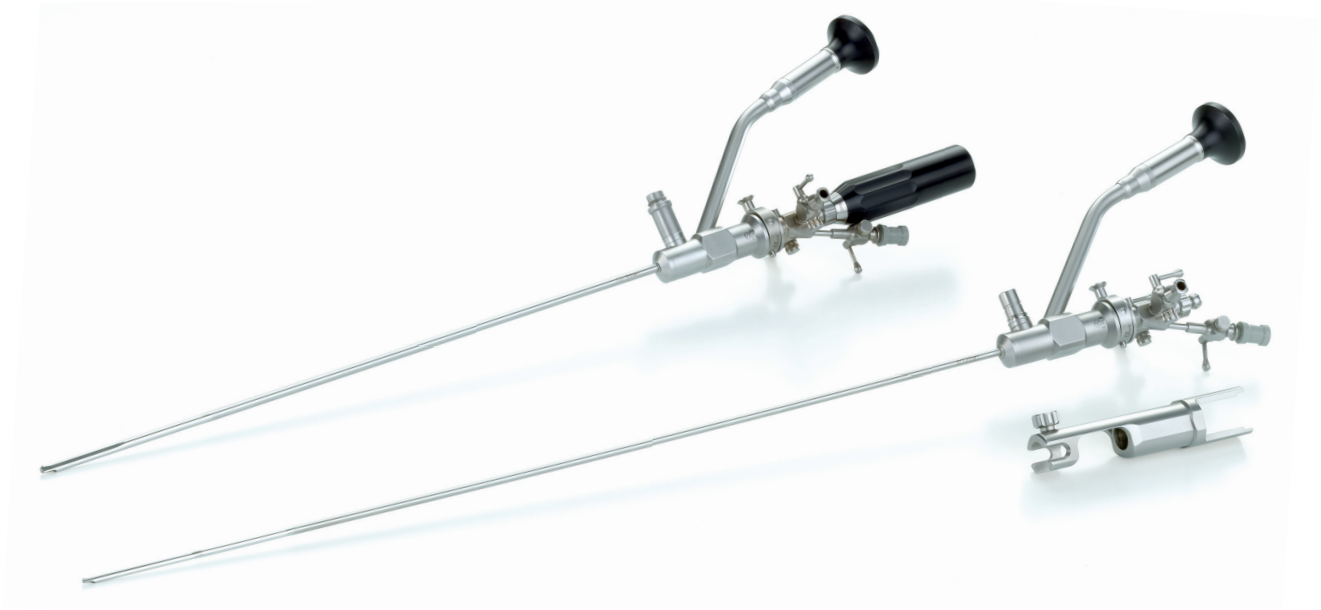

In order to insert the URS smoothly into the natural orifice of the urethra, the bladder and further into the ureter, the distal tip was realized in a special, atraumatic form. There is no need of camera re-focusing and no loss of picture while the URS is bended due to the ocular construction. Viewing direction to instrument channel as well as picture quality during all stages of operation and in all areas is ideal. Safe and ergonomic handling is offered due to an easy locking mechanism of the instrument bridge.

- Outstanding performance (picture quality, insertion, irrigation)
- Safe and reliable (robust, laser welded, easy to clean)
- Durability extends life time, resulting in low cost of ownership
- Easy to use due to unique locking mechanism
- Smooth insertion into urethra and ureter by special distal tip
- Straight channel for rigid lithotripsy probes, laser fibers, stone baskets or other instruments

## OEM PRODUCTS – PUT YOUR STAMP ON OUR EXPERTISE



Logo/Labeling



SCHÖLLY FIBEROPTIC GMBH  
Robert-Bosch-Strasse 1-3  
79211 Denzlingen, Germany

Phone: +49 7666 908-0  
medical-endoscopy@schoelly.de  
www.schoelly.de/ME

# Uretero-Renoscopes (URS)



## TECHNICAL DATA

| WORKING Ø      | WORKING CHANNEL FOR INSTRUMENTS WITH A Ø OF | WORKING LENGTH | DIRECTION OF VIEW     | FIELD OF VIEW | ITEM NUMBER  |
|----------------|---------------------------------------------|----------------|-----------------------|---------------|--------------|
| Distal 7.5 Fr. | 1 x 5.0 Fr. or 2 x 3.0 Fr.                  | 310 mm         | 8° to working channel | 85°           | 41.0611a.wol |
| Distal 7.5 Fr. | 1 x 5.0 Fr. or 2 x 3.0 Fr.                  | 425 mm         | 8° to working channel | 85°           | 41.0609a.wol |
| Distal 6.0 Fr. | 1 x 4.0 Fr. or 2 x 2.4 Fr.                  | 310 mm         | 8° to working channel | 100°          | 41.0615a.wol |
| Distal 6.0 Fr. | 1 x 4.0 Fr. or 2 x 2.4 Fr.                  | 425 mm         | 8° to working channel | 100°          | 41.0614a.wol |

Scope of delivery: Uretero-renaloscope, detachable instrument bridge, one pack of sealing caps and one pack of three slot membranes

## ACCESSORIES


| DESCRIPTION                                  | ITEM NUMBER  |
|----------------------------------------------|--------------|
| EMS LC2 adapter; for all URS mentioned above | 41.0610a.lc2 |
| EMS LC1 adapter; for all URS mentioned above | 41.0610a     |
| Instrument bridge                            | 41.0617a     |
| Sealing cap, package of 10 pcs.              | 41.0600d     |
| Three slot membrane, package of 10 pcs.      | 41.0593d     |

## DETAILS

- Highly optimized distal tip design
- High irrigation flow
- Autoclavable
- Semi-flexible sheath: 15° bending in all directions at full lifetime
- With image bundle (17,000/30,000 pixels)



Distal tip

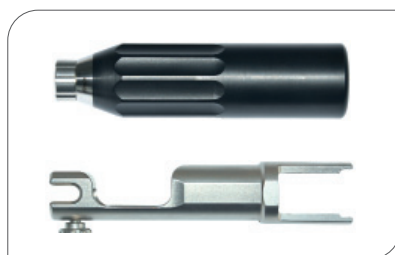

Adapters for lithotripter

Technical specifications are subject to change. The information is not legally binding. The contents are only for information about our products. Reprint, in whole or in part, is not permitted. Medical devices are only allowed to be marketed and used in countries where they have been approved for sales by regulatory authorities. For information on the available market access please contact us.

MEDICAL ENDOSCOPY  
MPA303-000-01-ME  
Version: B  
12-May-2021