

FLEXILUX II Laparoscopes

Ø 10 mm




The laparoscopes we have on offer are of high quality in optics and offer a superb contrast as well as an outstanding color reproduction. Equipped with high performance lenses they feature a full HD image quality.

As the illumination fibers are fastened strain free, they experience no tensile stress during the autoclaving process. Together with the welded outer tube, the system is sealed tight which results in a long life cycle for the optical components.


For laparoscopic applications in visceral surgery, urology and gynecology, the laparoscopes feature:

- Homogenous and sharp-edged illumination
- Optimized optical system for very bright images with rich contrast
- Strain free fastened illumination fibers
- Gold soldering of the plan glass which increases tightness and durability
- Low cost of ownership due to long life cycle

PRIVATE LABEL PRODUCTS - OUR EXPERTISE LABELED WITH YOUR NAME

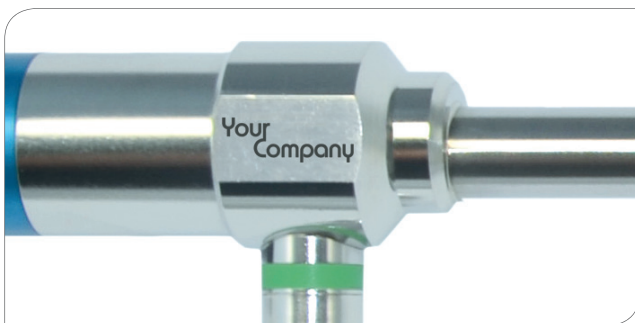

Logo/Labeling

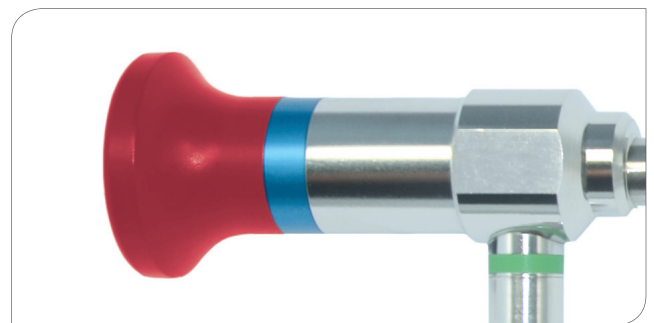

Color of eyepiece cap



SCHÖLLY FIBEROPTIC GMBH
Robert-Bosch-Strasse 1-3
79211 Denzlingen, Germany

Phone: +49 7666 908-0
Fax: +49 7666 908-380

www.schoelly.de/ME
info@schoelly.de

FLEXILUX II Laparoscopes

Ø 10 mm



TECHNICAL DATA

WORKING Ø	WORKING LENGTH	DIRECTION OF VIEW	FIELD OF VIEW	ITEM NUMBER
10 mm	344 mm	0°	wide angle	11.0031a
10 mm	344 mm	30°	wide angle	11.0043a
10 mm	344 mm	45°	wide angle	11.0041a

DETAILS

- State-of-the-art tip segment
- Quality feature: HDScope
- Autoclavable



Direction of view: 0°



Direction of view: 30°



Direction of view: 45°

Technical specifications are subject to change. The information is not legally binding. The contents are only for information about our products. Reprint, in whole or in part, is not permitted. Medical devices are only allowed to be marketed and used in countries where they have been approved for sales by regulatory authorities. For information on the available market access please contact us.

MEDICAL ENDOSCOPY
MPA301-000-01-ME
Version: 0
2018-Feb-06