

FLEXIVISION® 10

Controller Full HD



- Outstanding performance
 - Edge enhancement, automatic gain control (AGC) and dynamic contrast function for best in class performance
 - Latest sensor technology offers: excellent image quality at competitive pricing, no pixel errors due to automatic pixel defect compensation
- Multi-connectivity
 - Camera head full HD zoom
 - Camera head full HD fix focal length
 - Camera head full HD 90° angled
- Time saving and useful camera presets
 - Optimized settings for laparoscopy, arthroscopy, gynecology, urology, ENT, fiberscopes, neuroendoscopy
- SCHÖLLY MIS-Bus interface for direct light source control via camera head

PRIVATE LABEL PRODUCTS - OUR EXPERTISE LABELED WITH YOUR NAME



Logo/Labeling



Front frame color



Button color



PRIVATE LABEL

SCHÖLLY FIBEROPTIC GMBH
Robert-Bosch-Strasse 1-3
79211 Denzlingen, Germany

Phone: +49 7666 908-0
Fax: +49 7666 908-380

www.schoelly.de/ME
info@schoelly.de

FLEXIVISION® 10

Controller Full HD



TECHNICAL DATA

DESCRIPTION	ITEM NUMBER
Controller Full HD	95-3910
TECHNICAL DATA	
Operation control panel	On-Off, menu, white balance
Video outputs	High definition: 2 x DVI (1080p), 2 x HD-SDI (1080i)
Remote outputs	2 x 3.5 mm phone jack connectors (to connect recording/printer devices)
Remote inputs	2 x 3.5 mm phone jack connectors (to connect a foot control switch)
Connectivity	SCHÖLLY MIS-Bus interface
Dimensions/Weight	295 x 100 x 355 mm (W x H x D) / 4.39 kg
SYSTEM	
Standard compliance	IEC 60601-1 2 nd and 3 rd
Directive compliance	CE-mark
Power connection	100 - 240 V AC, 50/60 Hz
Power consumption	0.8 - 0.4 A

ACCESSORIES

DESCRIPTION	ITEM NUMBER
USB keyboard, German layout	PS00510
USB keyboard, US layout	PS00511
Foot control switch 1 pole, momentary mono jack plug with IPX7	6226-ACAA-BZZZ-002
Reprocessing tray for camera heads	95-4860


PLEASE SEE ADDITIONAL DATA SHEETS

FlexiVision Camera Head Full HD Zoom | FlexiVision Camera Head Full HD Fix Focal Length | FlexiVision Camera Head Full HD 90° Angled | FlexiLux 300 LED Light Source | FlexiLux 200 LED Light Source | FlexiVision 10 Controller Full HD | FlexiVision 20 Controller Full HD with USB Image/Video Capture | FlexiVision 20 Controller Full HD with USB Image Capture

DETAILS



Multi-connectivity

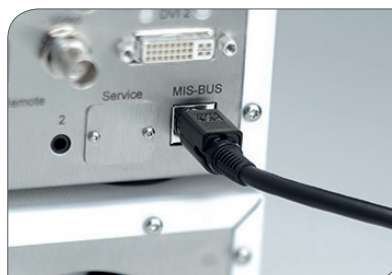

MIS-Bus



Accessories

Technical specifications are subject to change. The information is not legally binding. The contents are only for information about our products. Reprint, in whole or in part, is not permitted. Medical devices are only allowed to be marketed and used in countries where they have been approved for sales by regulatory authorities. For information on the available market access please contact us.

MEDICAL ENDOSCOPY
MPA310-000-01-ME
Version: A
2019-Apr-11