

FLEXILUX 4 Laparoscopes

Ø 10 mm



The new laparoscopes “FlexiLux 4” are optimized for the use with 4K UHD camera systems. Assembly, alignment and adjustment are carried out on a special qualified production line.

- High depth of field: A large depth of field range ensures sharp imaging even with changing working distances. Therefore refocusing during surgery is usually not necessary.
- Sharp images to the edge: The resolution has been improved especially in the edge areas resulting in a homogeneously sharp image in ultra-high resolution.
- High-contrast image: The optimized objective reduces image artifacts and reflections; thus maximizes the contrast.

Furthermore the laparoscopes feature:

- Homogenous illumination in the entire image area
- Strain free fastened illumination fibers
- Gold soldering of the plan glass which increases tightness and durability
- Low cost of ownership due to long life cycle

PRIVATE LABEL PRODUCTS - OUR EXPERTISE LABELED WITH YOUR NAME




Logo/Labeling



Color of serial ring



Shape of main body

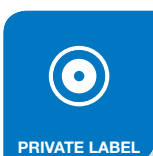

SCHÖLLY FIBEROPTIC GMBH
Robert-Bosch-Strasse 1-3
79211 Denzlingen, Germany

Phone: +49 7666 908-0
medical-endoscopy@schoelly.de
www.schoelly.de/ME

FLEXILUX 4 Laparoscopes

Ø 10 mm



TECHNICAL DATA

WORKING Ø	WORKING LENGTH	DIRECTION OF VIEW	FIELD OF VIEW	ITEM NUMBER
10 mm	344 mm	0°	wide angle	11-0231a
10 mm	344 mm	30°	wide angle	11-0243a

PLEASE SEE ADDITIONAL DATA SHEETS

FlexiVision 40 Controller 4K UHD Image Video Capture | FlexiVision Camera Head 4K UHD Zoom | FlexiLux 300 LED Light Source

DETAILS

- State-of-the-art tip segment
- Quality feature: UHDScope
- Autoclavable



Direction of view: 0°



Direction of view: 30°

A perfect match for perfect colors

BT2020 & CRI >90

Large camera color space BT2020 High CRI Ra >90

FlexiVision 40 camera with FlexiLux 300 LED light source

Technical specifications are subject to change. The information is not legally binding. The contents are only for information about our products. Reprint, in whole or in part, is not permitted. Medical devices are only allowed to be marketed and used in countries where they have been approved for sales by regulatory authorities. For information on the available market access please contact us.

MEDICAL ENDOSCOPY
MPA325-000-01-ME
Version: 0
2019-Oct-15