

Mini Endoscopes

Semi-flexible




Suitable for several applications in the surgical specialties:

- Arthroscopy
- Sinuscopy
- Laparoscopy

The semi-flexible mini endoscopes perfectly cover those fields where accessibility is difficult and space limited. Best examples are small joint arthroscopy or the access of blocked sinuses. In laparoscopy the mini endoscopes are used for special applications. E.g. for a view into the abdomen right from the first puncture the mini endoscopes can be used with a special optical veress needle. The semi-flexible mini endoscopes with high resolution image bundles are less prone to break due to their special construction, thereby considerably reducing repair expenses.

In combination with a FlexiVision camera system and the video algorithm "grid removal", the honeycomb structure produced by image bundle endoscopes can be visibly reduced, therefore further enhancing the image quality.


PRIVATE LABEL PRODUCTS - OUR EXPERTISE LABELED WITH YOUR NAME

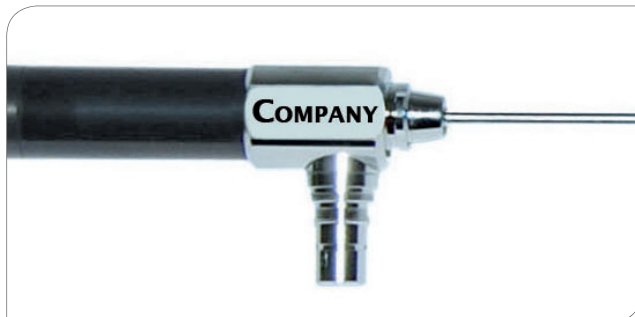

Logo/Labeling



Color of eyepiece cap

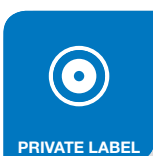

SCHÖLLY FIBEROPTIC GMBH
Robert-Bosch-Strasse 1-3
79211 Denzlingen, Germany

Phone: +49 7666 908-0
Fax: +49 7666 908-380

www.schoelly.de/ME
info@schoelly.de

Mini Endoscopes

Semi-flexible



TECHNICAL DATA


TYPE	ENDOSCOPE Ø	WORKING LENGTH	DIRECTION OF VIEW	IMAGE BUNDLE	ITEM NUMBER
Laparoscope	2 mm	258 mm	0°	50 k	11.0480a
Laparoscope	2 mm	258 mm	30°	50 k	11.0475a
Arthroscope	2 mm	110 mm	0°	50 k	21.0622a
Arthroscope	2 mm	110 mm	30°	50 k	21.0623a
Sinuscope	2 mm	110 mm	0°	50 k	31.0813a
Sinuscope	2 mm	110 mm	30°	50 k	31.0814a

ACCESSORIES


DESCRIPTION	ITEM NUMBER
Protection tube for endoscopes ø 4 mm and smaller	ye.0021a

DETAILS

- Autoclavable
- Semi-flexible sheath: 15° bending in all directions at full lifetime

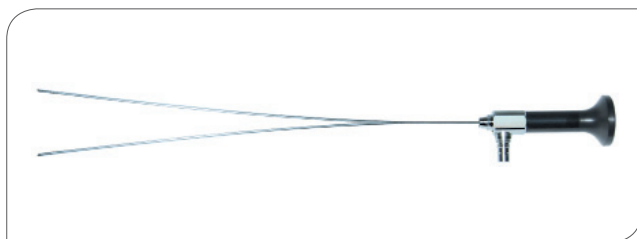

Semi-flexible

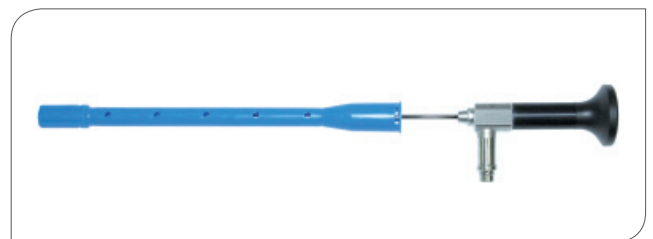

Protection tube

Technical specifications are subject to change. The information is not legally binding. The contents are only for information about our products. Reprint, in whole or in part, is not permitted. Medical devices are only allowed to be marketed and used in countries where they have been approved for sales by regulatory authorities. For information on the available market access please contact us.

MEDICAL ENDOSCOPY
MPA300-000-01-ME
Version: 0
2018-Feb-06