

# FLEXILUX II Arthroscopes

Ø 2.0 mm, 2.7 mm and 4.0 mm



The arthroscopes we have on offer are of high quality in optics and offer a superb contrast as well as an outstanding color reproduction. The optical system of the 2.7 mm and 4.0 mm scopes is calculated and aligned with regard to the best possible resolution. This ensures that it can record the HD resolution of the endoscopy camera.

The 2.0 mm arthroscopes are built with an image bundle. Due to the small diameter, they are especially suited for small joint arthroscopy.

The arthroscopes in 2.7 mm and 4.0 mm are mainly used for procedures on the knee, shoulder, jaw, wrist, ankle and elbow. At SCHÖLLY they feature:

- Optimized optical system for very bright images with rich contrast
- Low cost of ownership due to long life cycle
- Compatibility

In order to use sheath systems of other manufactures, SCHÖLLY offers arthroscopes which are mechanically compatible.

The outer dimensions and the optical parameters are matched to the use with the sheaths of the respective system/manufacturer.

## PRIVATE LABEL PRODUCTS - OUR EXPERTISE LABELED WITH YOUR NAME




Logo/Labeling



Color of serial ring



Shape of main body

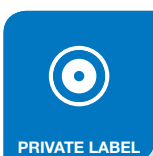

SCHÖLLY FIBEROPTIC GMBH  
Robert-Bosch-Strasse 1-3  
79211 Denzlingen, Germany

Phone: +49 7666 908-0  
medical-endoscopy@schoelly.de  
www.schoelly.de/ME

# FLEXILUX II Arthroscopes

Ø 2.0 mm, 2.7 mm and 4.0 mm




## TECHNICAL DATA

WORKING-Ø	WORKING LENGTH	DIRECTION OF VIEW	FIELD OF VIEW	TO SYSTEM	ITEM NUMBER
2.0 mm	110 mm	0°	wide angle	Storz	21.0622a
2.0 mm	110 mm	30°	wide angle	Storz	21.0623a
2.7 mm	187 mm	0°	standard	Storz	21.0617a
2.7 mm	187 mm	30°	standard	Storz	21.0619a
2.7 mm	187 mm	70°	standard	Storz	21.0621a
4.0 mm	175 mm	0°	wide angle	Storz	21.0003a
4.0 mm	175 mm	30°	wide angle	Storz	21.0008a
4.0 mm	175 mm	45°	wide angle	Storz	21.0101a
4.0 mm	175 mm	70°	wide angle	Storz	21.0011a
4.0 mm	170 mm	30°	wide angle	Wolf	21.0051a
4.0 mm	170 mm	70°	wide angle	Wolf	21.0053a
4.0 mm	142 mm	30°	wide angle	Stryker	21.0081a
4.0 mm	160 mm	30°	wide angle	Olympus	21.0066a
4.0 mm	157 mm	30°	wide angle	Dyonics	21.0061a

## DETAILS


- Quality feature for 2.7 mm and 4.0 mm arthroscopes: HDScope
- Semi-flexibility of 2.0 mm arthroscopes
- Autoclavable

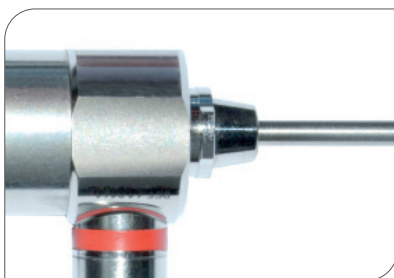

Connection system Storz



Connection system Wolf

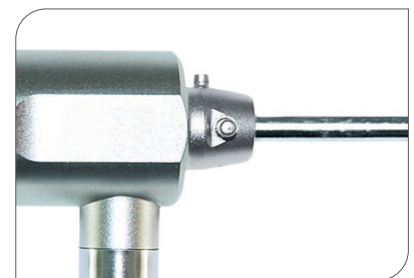

Connection system Olympus

Technical specifications are subject to change. The information is not legally binding. The contents are only for information about our products. Reprint, in whole or in part, is not permitted. Medical devices are only allowed to be marketed and used in countries where they have been approved for sales by regulatory authorities. For information on the available market access please contact us.

MEDICAL ENDOSCOPY  
MPA323-000-01-ME  
Version: 0  
2019-Aug-26